



Camelian Support



Catalogue No.	Proximal (F)	Distal (F)	Length (cm)	Inner Diameter Proximal / Distal (inch)	Tip	Hydrophilic Coating (cm)	Guidewire (inch)	Guide Catheter (inch)
---------------	--------------	------------	-------------	---	-----	--------------------------	------------------	-----------------------

CARNELIAN SUPPORT 14

MCP090A-E	2.6	1.8	90	0.021 / 0.017	Straight	60	0.014	0.035
MCP135A-E	2.6	1.8	135	0.021 / 0.017	Straight	100	0.014	0.035
MCP150A-E	2.6	1.8	150	0.021 / 0.017	Straight	120	0.014	0.035

CARNELIAN SUPPORT HD

MCPH090A-E	2.6	1.8	90	0.0205 / 0.017	Straight	60	0.014	0.035
MCPH135A-E	2.6	1.8	135	0.0205 / 0.017	Straight	100	0.014	0.035
MCPH150A-E	2.6	1.8	150	0.0205 / 0.017	Straight	120	0.014	0.035

CARNELIAN SUPPORT DP

MCPMS60A-E	2.6	1.8	60	0.0205 / 0.017	Straight	40	0.014	0.035
------------	-----	-----	----	----------------	----------	----	-------	-------

CARNELIAN SUPPORT BTA

MCPB150-E	2.6	1.6	150	0.021 / 0.0155	Straight	120	0.014	0.035
-----------	-----	-----	-----	----------------	----------	-----	-------	-------

CARNELIAN SUPPORT 18

MEP070-E	2.6	1.8	70	0.021 / 0.0195	Straight	40	0.018	0.035
MEP090-E	2.6	1.8	90	0.021 / 0.0195	Straight	60	0.018	0.035
MEP135-E	2.6	1.8	135	0.021 / 0.0195	Straight	100	0.018	0.035
MEP150-E	2.6	1.8	150	0.021 / 0.0195	Straight	120	0.018	0.035

FOR ANY CASE IN PERIPHERAL INTERVENTION

Camelian Support



1639

Tokai Support Catheter

Products listed may not be available in all markets. Please refer to the insert for complete instructions, contraindications, warnings and precautions.

TMP-1807284-03



Camelion

SUPPORT

Multistep Shaft with Optimal Balance

provides additional strength, pushability and crossability

Tungsten Braid and Platinum Marker Band

provide excellent visibility enabling precise positioning

Hydrophilic Coating

facilitates easy and smooth tracking

Extensive Selection

available for any case in peripheral intervention

The Proven Quality and Performance


the most accepted support catheter in Japan



Recommended Indications

	14	HD	DP	BTA	18
Antegrade approach	✓	✓			✓
Retrograde approach	✓	✓			✓
Trans-collateral approach	✓			✓	
SFA-CTO		✓			✓
Below the Knee				✓	
Direct approach			✓		

Stiffness

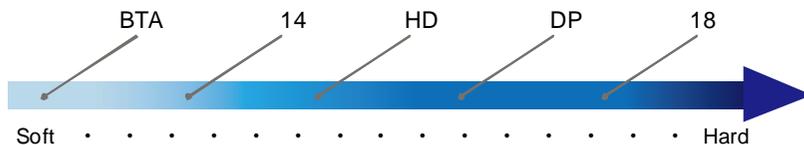


14 Support 14

Workhorse

1.8Fr/2.6Fr 90, 135, 150cm



- Well-balanced and Multi-purposed support catheter

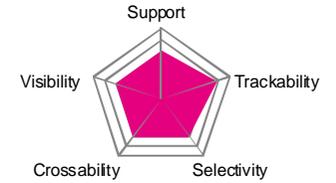


HD Support HD

Hard shaft

1.8Fr/2.6Fr 90, 135, 150cm



- Hard shaft provides additional back-up force

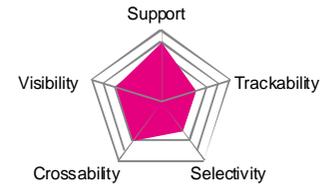


DP Support DP

Direct Puncture

1.8Fr/2.6Fr 60cm



- Specially designed for direct approach
- Short tapered tip provides easy insertion and prevents blood from leaking at puncture site

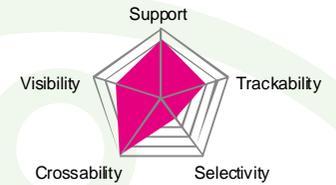


BTA Support BTA

Below The Ankle

1.6Fr/2.6Fr 150cm



- 300mm of 1.6Fr. flat shaft allows tracking through long and tortuous anatomy such as trans-collateral approach

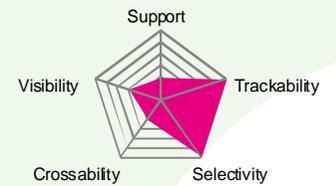

18 Support 18

0.018" compatible

1.8Fr/2.6Fr 70, 90, 135, 150cm



- Low profile tapered tip with additional strength of shaft facilitates easy crossing of challenging lesion.

